

# Segger Tutorial – SWIM (Segment-guided Water and Ion Modeling)


Last updated: Jan. 31, 2021 (Segger v2.5.4, Chimera Version 1.13)  
By: Greg Pintilie (gregdp@stanford.edu)

## Overview

- In this tutorial, we will use the **SWIM** dialog in **Segger** (from the Tools panel) to add Ions and Water molecules to a map of Apoferritin.
- Typically Ions and Water molecules can be resolved at  $\sim 2.5\text{\AA}$  resolution and higher for water and  $\sim 3\text{\AA}$  resolution and higher for ions. (These are very rough approximations and each situation can be different).


## Manual Modeling - Steps

1. Download or fetch the map EMD:22657, and the model, PDB:7k3v
2. Select and delete all existing ions and water molecules:
  - From the Chimera Select menu, select Residue, all non-standard.
  - From the Chimera Actions menu, select Atoms/Bonds, delete
3. Open the SWIM dialog from Segger Tools panel. It should look like this:

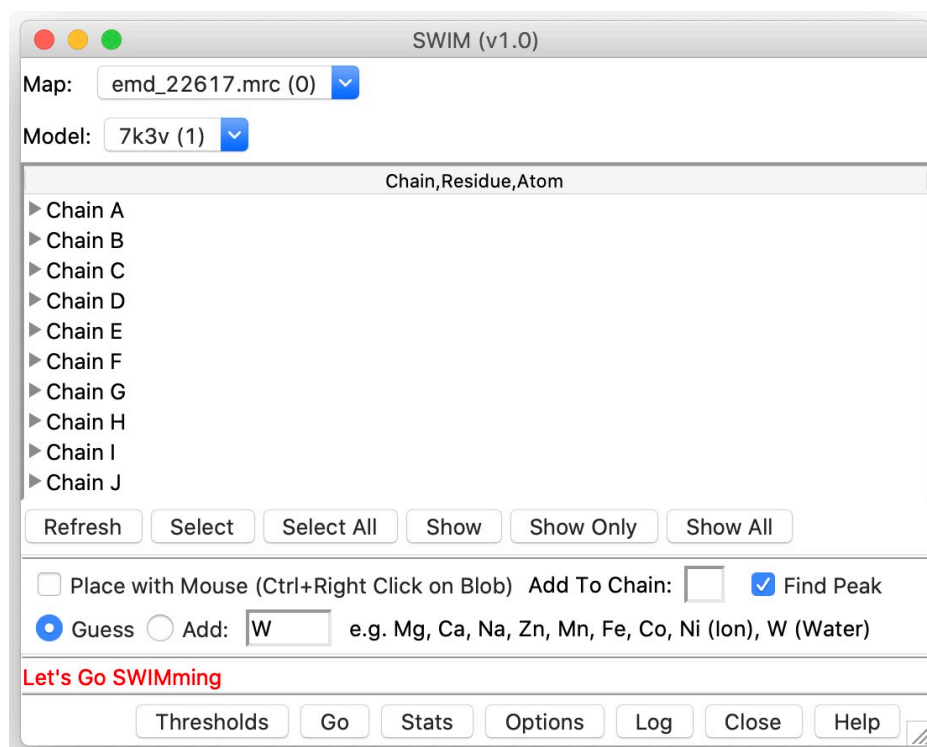


In the SWIM dialog, you can select a map and model to work on at the top. The list of chains are all the chains in the model. It is a tree list, so you can expand to see residues and atoms.


4. Select Chain A in the list, and press the **Select** button, then **Show Only** button under the list.

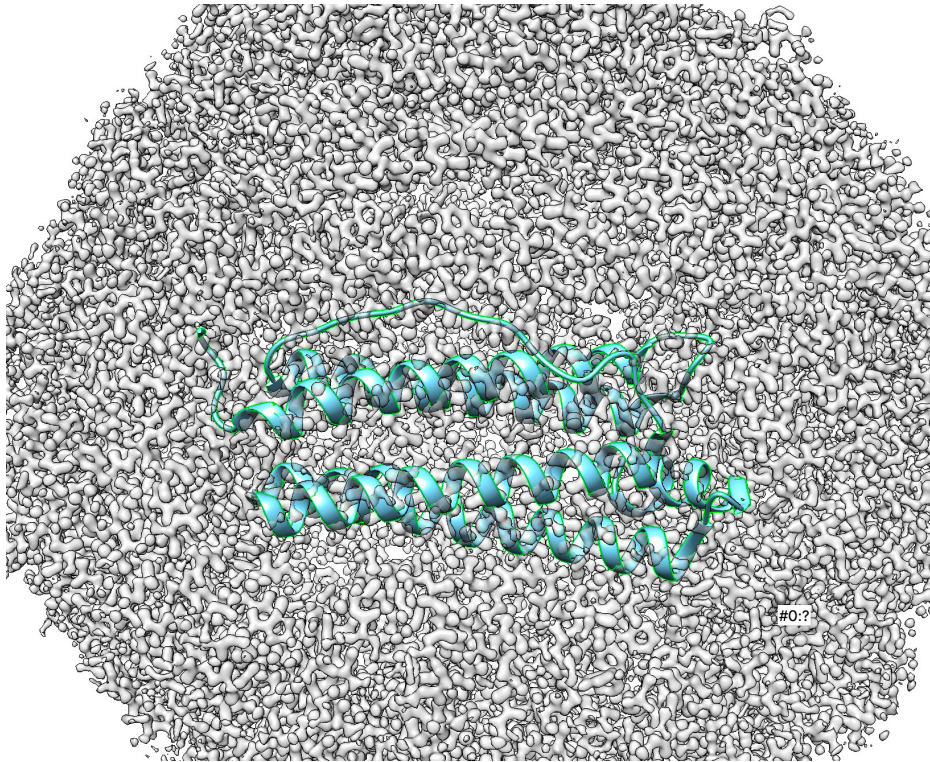


5. Press the **Zone** button in the SWIM dialog. This masks the full map with the selected atoms, with radius entered to the right of the Zone button (5Å by default). A new map will be created and shown (and the full map becomes hidden):

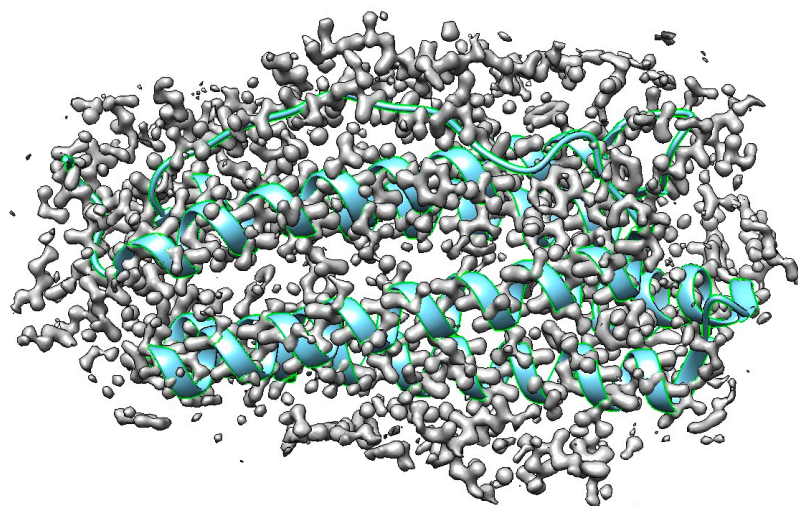

6. Now we can look up close for waters and ions. Show atoms using Actions menu on the Chimera window, Atoms/Bonds, show. Look around residue 27, and set the contour level of the zoned map to  $\sim 0.015$  in the Volume Viewer dialog.

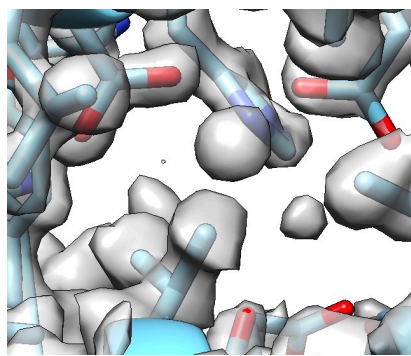

Before placing with Ctrl+Click

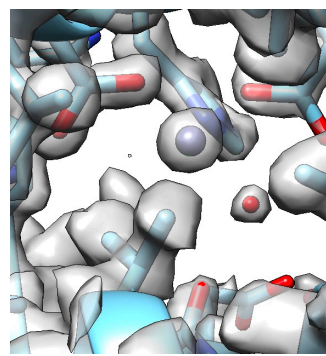

After

7. In the SWIM dialog, check '**Place with Mouse (Ctrl+Click)**', then Ctrl+Click on the blobs in near the residue. On the bigger blob, it will place an Ion, and on the smaller blob, it will place a water. This is if the '**Guess**' option is checked.

- To see the rationale, see the IDLE log (press Log at the bottom of the dialog to open it before placing the Ion/Water).
- The main criteria are the distances from the peak density in the blob to the nearest atoms, and the type of atom (polar or charged).
- If you press the **Options** button at the bottom of the SWIM dialog, the following panel is shown, where the distances can be adjusted:

Distance ranges (in Angstroms):

Ion distances: from  (Angstroms) to  (Angstroms)  
 Water distances: from  (Angstroms) to  (Angstroms)

8. When the '**Find Peak**' option is checked, the highest value in the 'blob', or the true peak of the blob is found (by uphill climb). If the option is not checked, the water/ion is placed on the highest value under where the Click is performed.

9. If '**Add: [ ]**' is checked, the atom entered in the text field is added instead of the '**Guess**' mechanism. Use Mg, Ca, Na, etc. for ions, and W for water.

10. '**Add to Chain [ ]**' may be changed to another letter or number. If left empty, the water/ion is placed as a residue in the same chain as the nearest atom.

11. When done adding manually, uncheck the '**Place with Mouse**' in the SWIM dialog. Note that Ctrl+Click normally is used to select atoms/maps, etc., so this brings back this standard functionality.

12. Now let's try placing more ions/waters **automatically** using the **Go** button.

- One consideration is what contour level to use in the map (above 0.015 was suggested). With the full map selected in the SWIM dialog, press the **Thresholds** button at the bottom. The contour levels are shown in the status line (in red) in the SWIM dialog and in IDLE, as follows:

```
1-sigma=[0.0068]
2-sigma:[0.0136]
3-sigma:[0.0204]
-- in emd_22617.mrc
```

- 1-sigma represents 1 standard deviation above the average of all values in the map, 2-sigma 2 standard deviations, and 3-sigma 3 standard deviations.
- A sensible value may be 2-sigma above the mean, to ensure that waters/ions are placed on peaks that are significantly above the background noise.
- A second consideration is whether the peak is well-resolved, which typically means its profile should be closer to a Gaussian function. To avoid placing waters/ions on less resolved peaks, the following option can also be checked in the **Options** panel in the SWIM dialog:

☐ Put water/ion only when Q-score >  sigma:


- The higher the Q-score, the more resolved the peak has to be, or closer to a Gaussian with a sigma (~width) of 0.4Å. Here we use sigma of 0.4 for maps up to 1.2Å maps; for resolutions up to 1.5Å, sigma of 0.6Å should be used. This map is ~1.3Å so we use 0.4.

13. In the Volume Viewer dialog, enter the contour level for the zone map (with the **\_Z5\_** in the name) as 0.0136 (2-sigma).


- When zoning the map for water/ion placement, make sure to use a zone radius of at least 5, so as waters/ions can be found at least 3.3Å away plus the drop-off on the other side of the peak.



14. In the Segger dialog, select the emd22657\_Z5\_map in the 'Map:' field, enter 0 for the number of steps in either 'Group by Connectivity' or 'Group by Smoothing' (we want to keep the original watershed regions without grouping), and press the **Segment** button. You should see something like this (you can select all the regions and press the **Tr.** button in the Segger dialog to make them transparent):

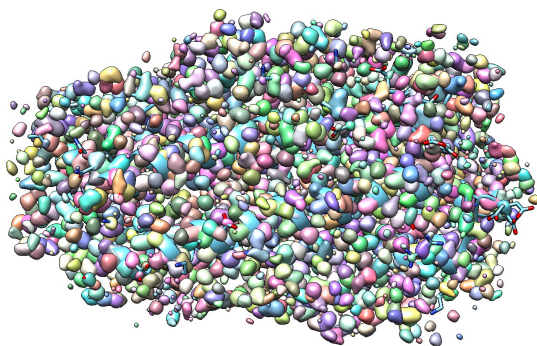


Regions after segmenting

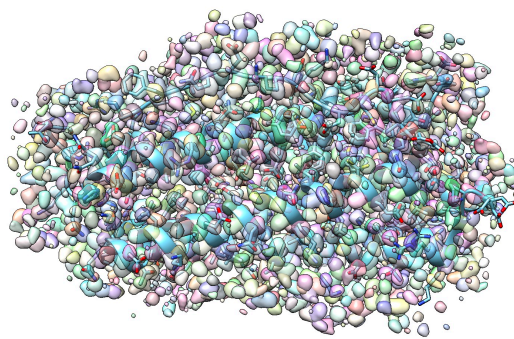


Regions after making them transparent with the Tr. button

15. Back in the SWIM dialog, select Chain A in the list, and press the **Select** button. The automatic placement only looks for peaks near selected atoms. Then press the **Go** button. You will see a progress in the SWIM dialog status and the Chimera window status. After about a minute, you should see the results, 198 waters and 31 ions. The process makes only the regions used for placement visible, and hides the rest.

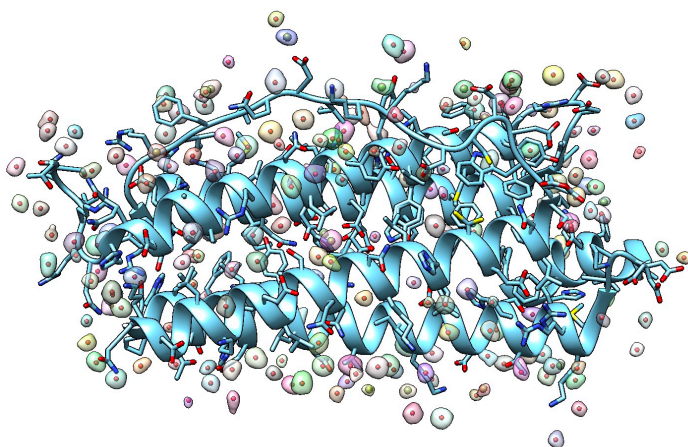


All waters and ions

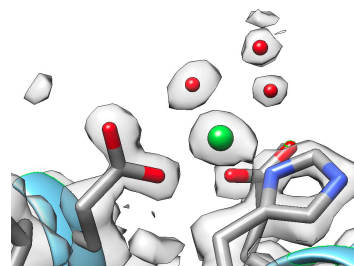


A close-up on Glu27

16. Finally, press the **Stats** button at the bottom of the SWIM dialog. This will collect some statistics on the distances between the waters/ions and the nearest atoms in the model, and also on Q-score distributions. The output is written to the IDLE window

(make sure to open the **Log** before pressing Stats). This output can be copied and pasted to Excel, for example, to create plots like this:

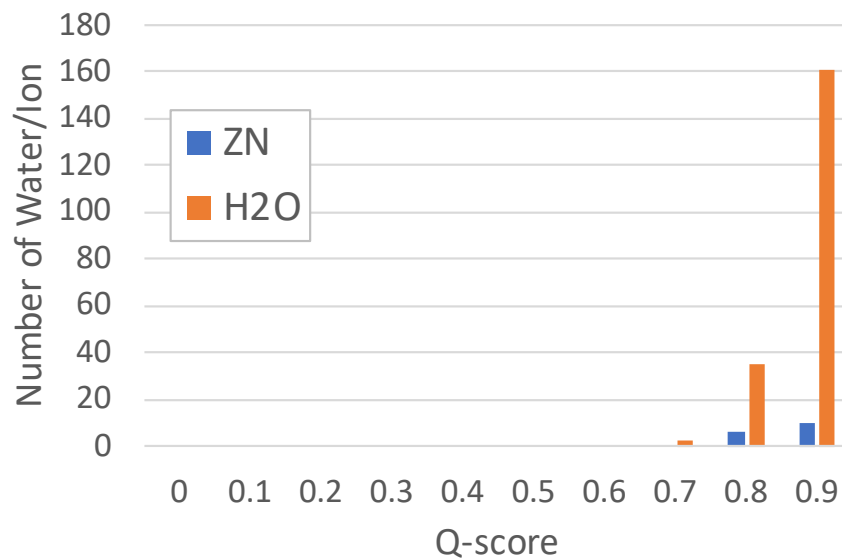

- The above shows the distribution of Q-scores for placed ZN and H2O atoms. Note that the automatic placement defaults to using Zn for atoms. This default can be changed by placing another ion (e.g. Ca) in the “Add: [ ]” field in the SWIM dialog.

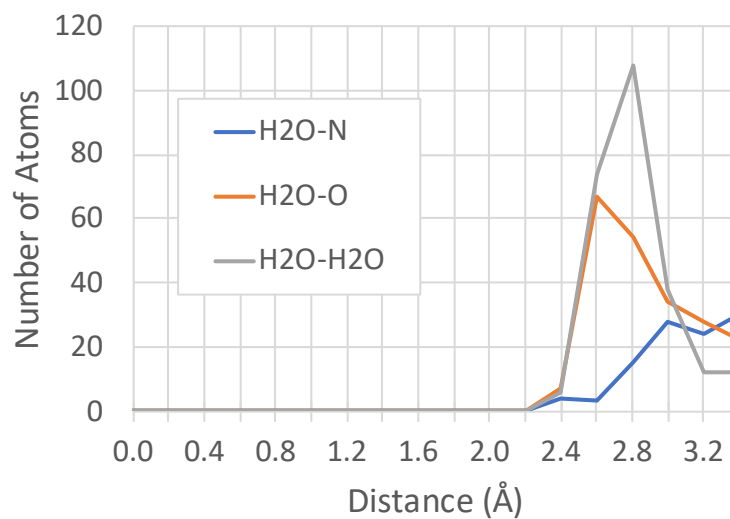

- The above shows a radial-distance plot for placed water atoms to nearby atoms of type N or O. Peaks for H2O-H2O and H2O-O can be seen at  $\sim 2.8 \text{ \AA}$ , which is the typical distance for hydrogen bonds.