

Segger Tutorial – Segger Dialog


Last updated: Jan. 31, 2021 (Segger v2.5.4, Chimera Version 1.13)

By: Greg Pintilie (gregdp@stanford.edu)


Opening a density map

As an example, the density map of [GroEL](#) at 10.0Å resolution (EMD:1080), from the [EMDB](#), will be used.

- It can be opened in Chimera by selecting File->Open....
- Instead of downloading it, you can also more simply choose "Fetch by ID ..." from the File menu, select "EMDB" and use the ID 1080. Either way, once open, the density map will appear in the main Chimera window:

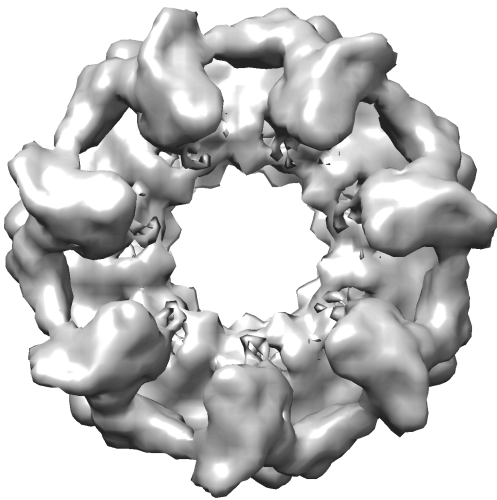


EMDB: 1080 as it appears once opened

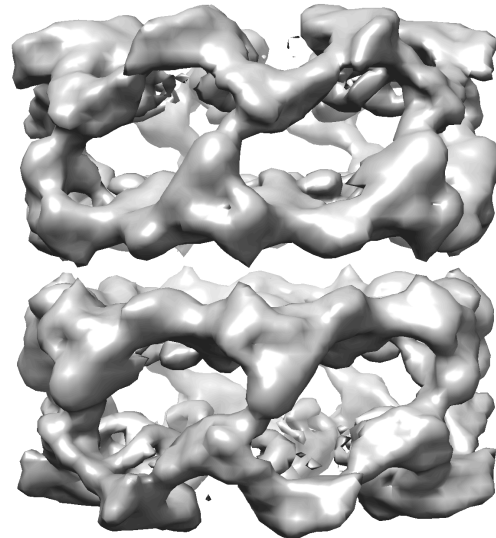

EMDB:1080 rotated around the x axis by ~90°

The Volume Viewer dialog will also appear if it wasn't shown already:

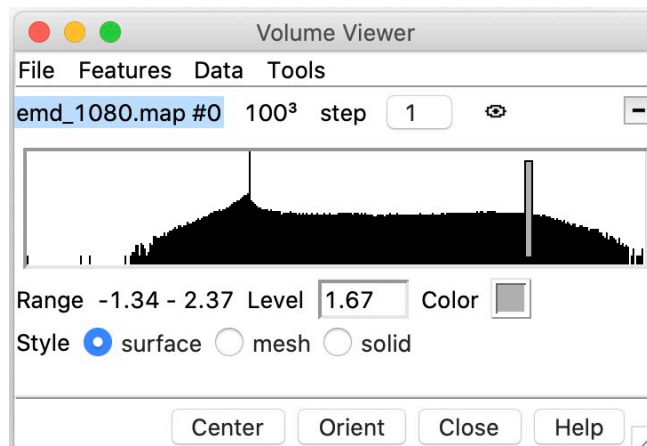

The Volume Viewer dialog can also be opened from Tools / Volume Data / Volume Viewer. It shows, amongst other information, a histogram of density values in the map. The gray bar can be dragged left or right to change the *threshold*, which is the density value used to generate the iso-surface shown in the main window. At higher thresholds, the inner and denser regions are seen, and at lower thresholds, the outer surface of the molecular complex is seen.

Opening the Segger dialog

The **Segger** dialog can be opened by selecting Tools / Volume Data / Segger. It looks like this:




Segmenting the map

To segment the map, first make sure the map to be segmented is selected in the field to the right of “Map:” in the Segger dialog. This field is a drop-down menu, showing all the density maps currently open.

- *Note* that if you open the Segger dialog after opening one or more density maps, an open map will already be selected by default. To select another map, click on the field, and then select it from the drop-down menu that appears.

When segmenting a map, a density threshold is applied, and only density values above this threshold will be included in the segmentation. The threshold chosen for the density map in the Volume Viewer interface is used; this is the same threshold (or contour level) used to draw the surface you see in the main Chimera window.

Once you've adjusted the threshold, to for example **1.12** as in this example, click the **Segment** button. After a few seconds, you will see the segmentation results in the main Chimera window, replacing the map previously shown there:

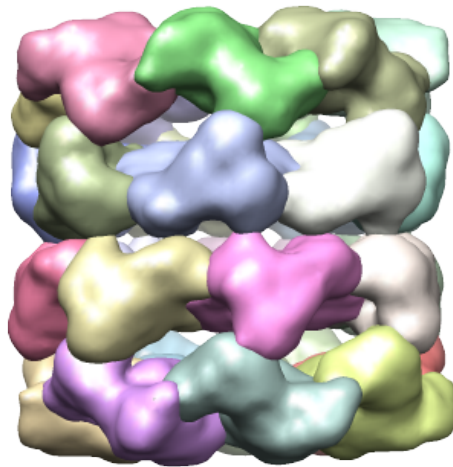

In the above segmentation, groups of two regions correspond to each protein. To see how to get a segmentation in which each region matches a single protein, read on...

(The **Center** button to the right of the map centers the camera rotation on the selected map).

Results of segmentation

Segmentation produces a number of regions, with each region containing only adjacent (26-connected) voxels. Each region is shown as a smooth surface surrounding the voxels it contains.

All the segmentation information is kept in a separate model in Chimera, which for the example above will have the name `emd_1080.seg`. You will see this model appear in the Model Panel dialog after pressing the "Segment" button.

After segmenting, also notice three changes in the Segger dialog:


- The name of the segmentation output file will appear in the field to the right of *Current segmentation*. In this example, it's `emd_1080.seg`.
- To the right of that, it will say 28 regions. This is the number of regions produced.
- At the bottom of the Segger dialog, just above the Segment button, it will also say *428 watershed regions, smoothed 3 voxels to get 28 regions*. Look at this line for other messages as you segment.

Segmentation method


Segger uses the [watershed method](#). The regions produced by the watershed method are often small and numerous - an effect often referred to as *over-segmentation*. To reduce this effect, and produce larger and more meaningful regions, Segger uses [scale-space](#)

[filtering](#) to group the many, small, resulting regions, into larger, more meaningful regions.


The watershed regions produced for the GroEL map are shown below in first image on the left. To produce fewer and larger regions, the density map is smoothed and the regions are grouped based on the smoothed map. This is repeated for a number of steps. The step size specifies how much smoothing takes place. By default, 3 steps of size 1 are used, leading to a segmentation as shown in the image in the middle.

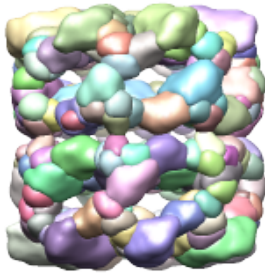

Watershed Regions

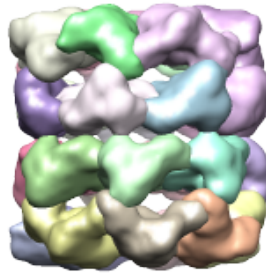

Regions after 3 steps of smoothing and grouping

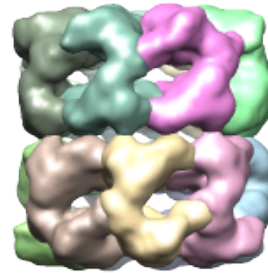


Regions after 6 steps of smoothing and grouping

There are several ways to direct this smoothing and grouping process:

- Continuing with the example so far, simply press the **Group** button. When no regions are selected, this button simply does more smoothing and grouping of regions. It does more smoothing steps automatically, until the number of regions changes. In this case, it takes 3 more steps, since only after this many more steps does the number of regions changes from 28 to 14.
- Click the **Options** button in the Segger dialog, enter 6 for the number of *Smoothing steps* (leaving 1 as step size), and then press the **Segment** button. Now the map is segmented again (which wasn't really necessary, but it does do a "clean" start in case you've been changing regions yourself already). Then 6 (this is the number you entered) smoothing steps are performed, and because more smoothing is performed, fewer regions are obtained - 14 regions. In this case, this is the same number of regions as proteins. The number 6 was arrived at after some experimentation, so do experiment with different numbers for your map.
- If we already knew there were 14 proteins in this density map, we can also enter a larger number of smoothing steps (e.g. 100), and 14 in the field to the right of *Stop grouping if number of regions <=*. With these options, smoothing steps (up to 100 of them) are taken, until the number of regions drops to 14. Of course it doesn't need to do all 100, it only does 6 as the other examples have shown.

Segmentation options

To change parameters for the segmentation process, press the **Options** button at the bottom of the Segger dialog to display the Options pane. The Segger dialog with the **Options** shown looks as follows:

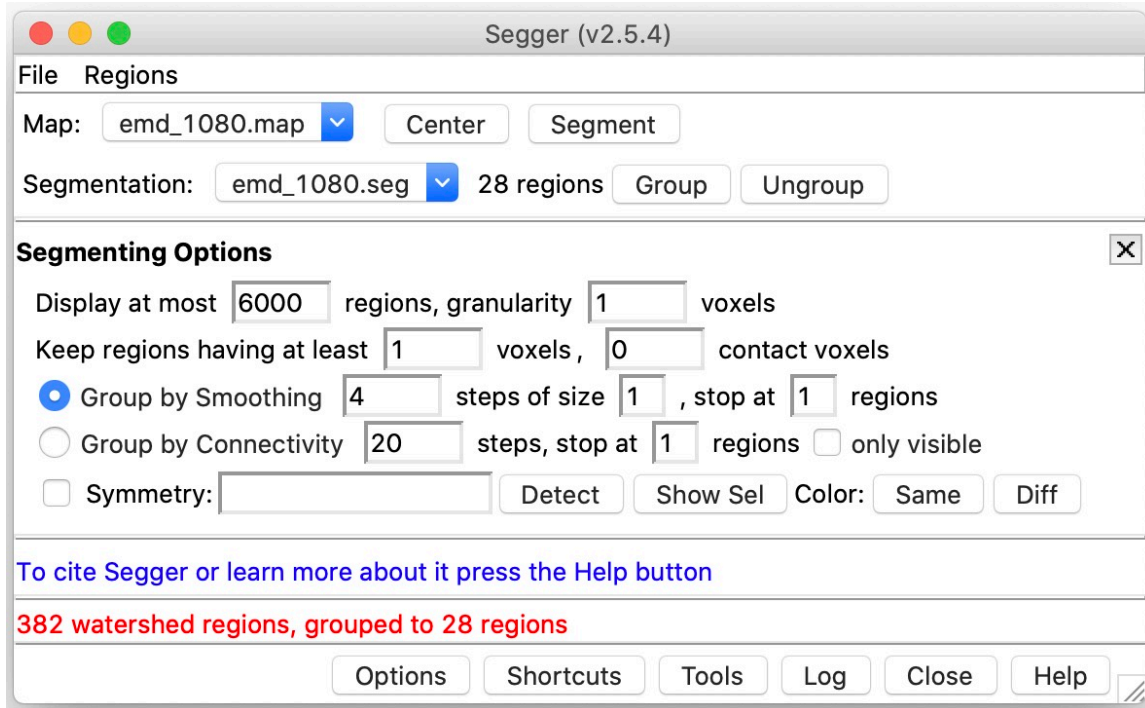

Display at most __ regions, granularity __ voxels

- The number of regions in some density maps can be quite large (tens of thousands). With this many regions, performance becomes very slow. Hence, the number entered in *Display at most ... regions* limits how many regions are shown. By default this number is 6000. The regions are still there in the segmentation, but they are just not shown.
 - If you segment a map and don't see regions where there should be regions, you can try increasing this number and then segmenting again.
 - Even when not all the regions are shown at first, they are still used when smoothing and grouping. So if you get more regions than you can view, you can try using more smoothing steps, which will produce fewer regions.
- Use larger granularity to speed up the display when there are many regions.

Keep regions having at least __ regions, __ contact voxels

- These options can be used to filter the resulting regions. Use higher # of voxels to remove small regions (typically noise). The granularity makes regions more coarse looking to speed up display when there are many regions present.

Group by Smoothing

- If this is selected, and 1 or more steps are used, then the watershed regions are grouped by smoothing the map in steps. If 0 steps are entered, then smoothing and grouping is not applied and the resulting regions are the watershed regions.
- You can specify a different number of smoothing steps and step size. If it looks like the regions obtained should be larger and there should be fewer of them, increase the number of smoothing steps.
 - Instead of changing this parameter, you can always:
 - Backtrack after segmenting – press the **Ungroup** button with no regions selected. This is the same as selecting one less number of steps and pressing segment again.
 - Group further by pressing the **Group** button with no regions selected. This is the same as entering one more than the number in the number of steps and pressing segment again.
- The step size parameter will also change the results, but not as significantly:
 - When smaller step sizes are used, less smoothing is performed at each step, and so the number of and the shape of each region will change more gradually; thus to get fewer regions, more steps will be required.
 - With larger step sizes, changes will be more drastic, but fewer steps will be required to get to a particular number of regions. Experiment with different numbers on your map... When the map contains a lot of noise, larger steps are recommended, since smoothing tends to reduce the bad effects of noise.
 - For ungrouping and regrouping purposes (discussed further below), use smaller step sizes and more steps, since this will give you a more gradual change in number of regions at each ungrouping step.
- The process also stops if it reaches the number of regions entered in stop at ____ regions.

Group by Connectivity

- When this is selected, watershed regions are grouped by looking at contacts between them. This tends to work better in higher resolution maps.
- Use fewer or more steps, or the Group and Ungroup buttons without regions selected.

Using Segger with your own map

Now that you have a better idea of how Segger works, and what the parameters do, go ahead and try Segger with your own map. Load the map in Chimera, select it in the field to the right of “Map:”, choose a threshold in the Volume Viewer dialog, and then press **Segment**. Adjust the segmenting parameters and try again, as necessary.

Selecting and working with regions

In the main Chimera window, where the segmented regions are shown, a region can be selected using the usual Chimera selection mechanism. For example, while holding the *Control/Ctrl* key, clicking on a region selects it (you will see a greenish border around it), or unselects it if it is already selected.

Multiple regions can also be selected; clicking on a region while holding *Shift+Control* adds it to the current selection if it isn't already in the selection, or removes it from the selection if it is in the selection already.


Other important operations are grouping and ungrouping regions.

- To group two or more regions, select them and then press the **Group** button at the bottom of the Segger dialog. The regions will now be combined, and shown using a single surface.
- To ungroup one or more regions, select it/them and then press the **Ungroup** button at the bottom of the Segger dialog. You can ungroup regions that you previously grouped, or which were grouped by smoothing and grouping! In the latter case, the effect is to 'undo' smoothing in select parts of the map, i.e. the parts that are covered by one or more regions...

Editing the segmentation

Grouping and ungrouping can be used to modify the resulting segmentation, to adjust for example a region that looks wrong. This process is illustrated below - first the orange region in the image on the left is ungrouped, producing smaller regions, which are then regrouped with nearby regions.


If you try to ungroup a region you obtained by smoothing and grouping (simply having pressed the Segment button with multiple smoothing steps), and it gets ungrouped into too many small regions, then you can try redoing the process, but using more steps, and steps that are smaller in size. This way, the grouping will be more gradual, and the intermediate regions sizes and numbers will drop more gradually.

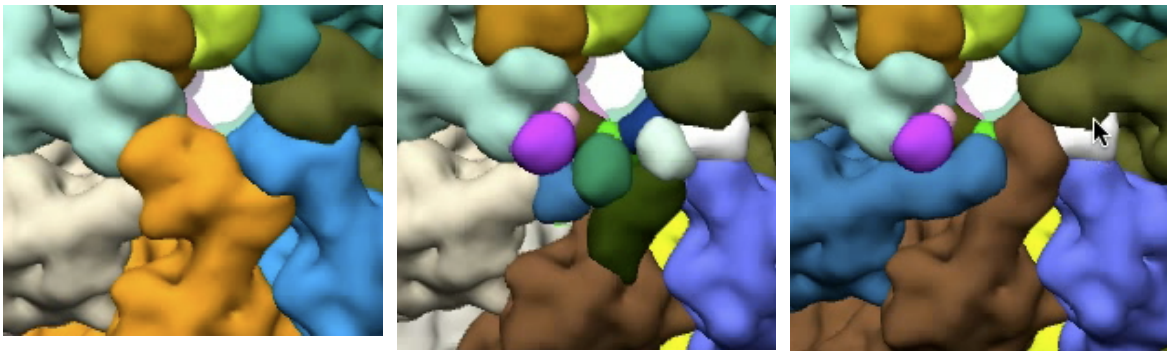

Shortcuts

Selected regions can be hidden, deleted, made transparent, etc., by selecting the one of the entries in the **Regions** menu which appears at the top of the Segger dialog. For convenience, some of these functions are accessible by pressing “Shortcuts” at the bottom of the Segger dialog:

Shortcuts

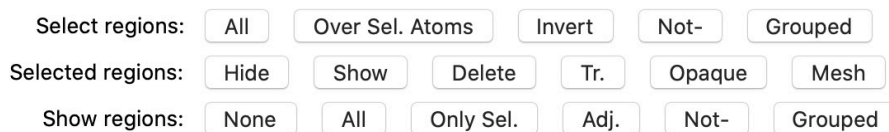

- **Select regions:**
 - Pressing one of these buttons will select regions in the current segmentation (the one selected at the top of the Segger dialog).
 - All: all regions
 - Over Sel. Atoms: regions that overlap selected atoms
 - Invert: all regions currently not selected
 - Not-: Regions that are not grouped
 - Grouped: Regions that are grouped
- **Selected regions:**
 - Pressing one of these buttons will perform an operation on the selected regions:
 - Hide: hide selected regions
 - Show: show selected regions
 - Delete: regions will be deleted from the segmentation
 - Tr.: selected regions become transparent
 - Opaque: selected regions become opaque
 - Mesh: display region as mesh instead of solid surface (use Tr. Or Opaque to change back to solid surface)
- **Show regions:**
 - Pressing one of these buttons will affect which regions are being shown
 - None: hides all regions
 - All: all regions, up to the number entered in Options for Display at most ____ regions. Note that regions are shown in order of size from largest to smallest, so the ____ largest regions are shown – if regions are interactively re-grouped, pressing this button is useful to reorder and show the new ____ largest.
 - Only Sel.: only the selected regions.
 - Adj.: regions that are adjacent to any selected regions.
 - Not-: regions that are not grouped
 - Grouped: only regions that are made up of 2 or more subregions.

Other Tools

Press the Tools button at the bottom to show some other tools:

- Extract: extract densities (or mask the map) using segmented regions
- SegFit – rigidly fit models based on segmented regions
 - Note this is the same as the dialog that shows up when selecting Tools -> Volume Data -> SegFit.
- rSeg – group regions based on radial distance from center
- iSeg – group regions based on icosahedrally-adjusted distance from center
- ProMod – combine multiple (flexible) fitting results into probabilistic models
- Movie – Make animations of models
- MapQ – Calculate Q-scores and visualize parts of a model
- SWIM – Segment-guided Water and Ion Modeling

Saving segmentation results

You can save and load a previously saved segmentation through the File menu at the top of the Segger dialog. Region colors are also saved along with the segmentation.